



Alcohol Pad

70% isopropyl alcohol


Alcohol pad is produced by acupuncture non-woven and 70% isopropyl alcohol which has slightly irritating smell. The product is flake type and packaged with aluminum foil bag. One pad saturated with 70% isopropyl alcohol. Single use only which is safe and continent. The product will not cause adverse stimulation to human skin and mucosa, will not cause allergic reactions and other damage, and will not precipitation toxic substances either.

Intended Use

Used to disinfect hard surface of non-invasive medical device.

Specifications

- Item: Alcohol Pad
- Material: Acupuncture Non-woven 130g/ m²
- Content: 70% Isopropyl Alcohol
- Function: Disinfection
- Pack Size: 50*50mm
- Pad Size: 30*30mm
- Feature: *Single layer, no folding
*Aluminum film bag 73g/ m²

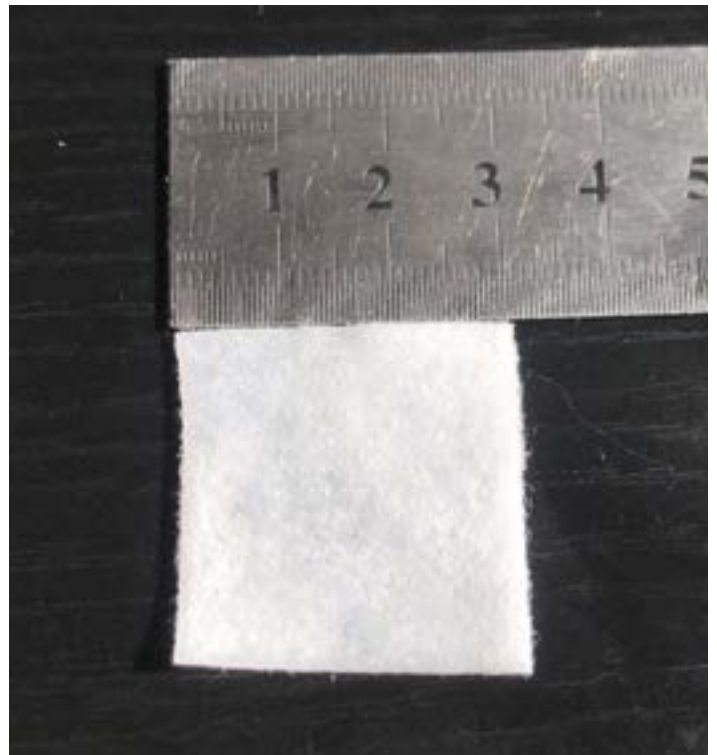



Directions

- One pad saturated with 70% Isopropyl Alcohol.
- Rub skin briskly in a circular motion from injection and injury site outward.
- For External single use only.
- Discard after single use.

Caution

- Not used to disinfect invasive medical device.
- Do not use when the package is unpack or damaged.
- In the event of skin irritation or rash,discontinuance
- Avoid contact with eyes,keep away from children.

Storage and Transportation

1. The device shall be stored at a relative humidity of less than 80%, and the storage environment shall be kept non-corrosive gases, dry, avoiding sunlight, well-ventilated conditions, and clean.
2. The device shall be taken care of during transportation and handling.

Packing Details

Package:	100pcs/box,100boxes/ctn
Carton size:	53.5*32*28.5cm
G.W/N.W:	13/12kg

