



## Three-way Stopcock


## Disposable

Three-Way Stopcock is mainly used for regulating the flowing rate of infusing fluids when clinical infusion.

### Indications and uses

The Disposable Stopcocks are intended for single use of intravenous infusion therapy, for changing the infusion direction(s), extending the infusion pipeline, and (or) pressure monitoring.

The clinical use time for the stopcock extension is not more than 7 days, and the clinical use time for other models of the product specifications is not more than 24 hours.

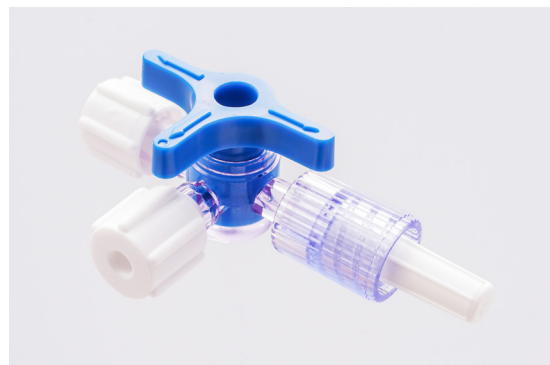

### Specifications

· Item	Three-Way Stopcock
· Material	Medical-grade polycarbonate material with good lipid resistance
· Options available	Two/Three Way/With or without extension tube
· Extension Tube	3.6mm(OD)*2.2mm(ID)*(5-300cm)
· Sterilization	EO sterilization
· Risk-related classification	Class II a, according to MDD Annex IX, rule 2
· Certificate	CE/ISO13485/FSC
· Shelf Life	3 years
· OEM	Available
· Sample	Available



## Contraindications


- Single use only, do NOT reuse. Reuse can lead to serious inflammation or infection even death.
- Dispose of the device as a contaminated medical device and in accordance with hospital protocols.

## Warning and Caution

No absolute contraindication.


## Main Parts of the device

### (1) Three-way stopcock

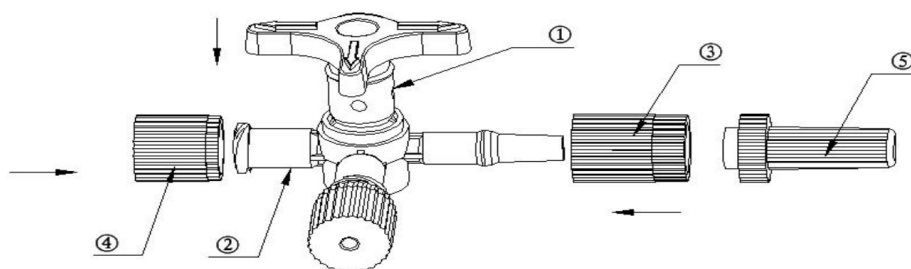

1 Core 2 Base 3 Lock 4 Male Protective cap 5 Female Protective cap

### (2) Three-way stopcock with extension tube

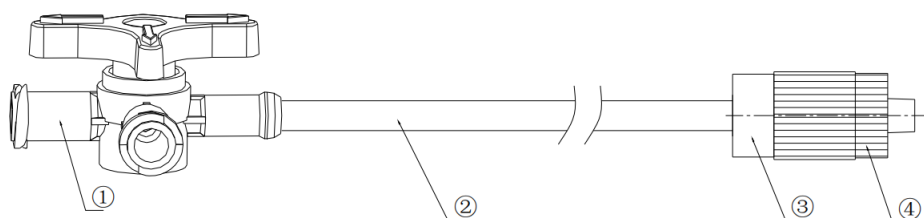

1 Stopcock 2 Extension tube 3 Joint 4 Lock





## Material list

Parts	Materials	Blood-contacting
Core	PE	Yes
Male Protective Cap	PP	No
Female Protective Cap	PP	No
Lock	PC	Yes
Base	PC	Yes
Tube	TPU/PVC	Yes

## Packing Details

The Disposable Stopcocks is packaged in primary packaging, secondary packaging, and tertiary packaging.

Primary packaging: Composite film blister heat sealed with TYVEK paper.

Secondary packaging: Single layer corrugated paper.

Tertiary packaging: Double walled corrugated board.

### (1) Three-way stopcock

Package: 50pcs/box, 20boxes/ctn

Carton size: 440\*340\*530mm

G.W/N.W: 8.6/8kg



### (2) Three-way stopcock with extension tube

Package: 100pcs/box, 4boxes/ctn

Carton size: 610\*350\*580mm

G.W/N.W: Calculated by actual extension tube ( 5-50cm )

